


A.5660

A.5650



Legatore Emorroidario Monouso LEM

LEM Disposable Haemorrhoid Banding Ligator

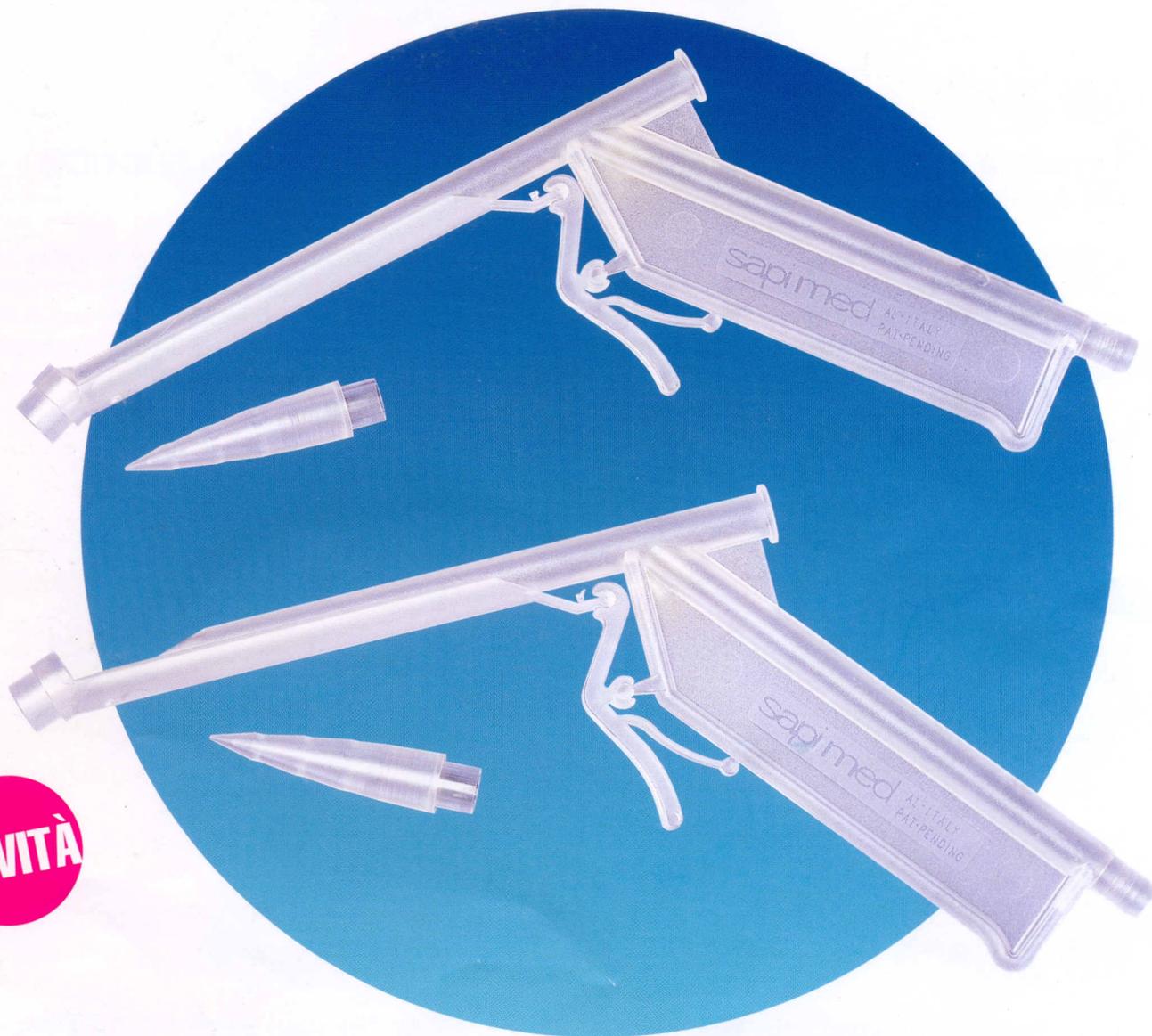

NOVITÀ

Il legatore emorroidario monouso "LEM" è un dispositivo SAPI MED brevettato ed unico nel suo genere.

È realizzato in due versioni: ad "ASPIRAZIONE" e a "TRAZIONE". Entrambi sono dotati di ampio cilindro (Ø 10mm) per il contenimento del nodulo. Dispone di un'impugnatura ergonomica ed il grilletto di comando espulsione anello elastico è provvisto di balestra di richiamo. Una calcolata angolatura permette una buona visione nell'intervento. Il cono di carico anellini è stato studiato per facilitare al massimo l'operazione. Il dispositivo è costruito in materiale plastico trasparente antiurto e marchiato CE secondo la direttiva 93/42.

The Sapimed 'LEM' Disposable Haemorrhoid Banding Ligator is a unique, patented instrument. It is manufactured from plastic, designed for single-patient use, and available in two variants: one for use with suction, and the other for use with forceps. Both variants are supplied with a separable cylindrical cone (Ø 10 mm) to assist in loading the bands, an angled distal cylinder to contain the nodule, and an ergonomically-designed pistol handgrip with trigger to release the banding ring. The LEM's angled design allows the operator a clear view of the procedure. The LEM banding ligator is designed for use with Sapimed Self Light operative proctoscopes.

